



U S E R M A N U A L



CONTENTS

CONGRATULATIONS ON YOUR VELA BLUES 100	3
1.0. OPERATION	4
1.2. BEFORE OPERATION	4
1.3. OPERATION	5
2.0. SAFETY POINTS	8
3.0. SET-UP	9
3.1. BATTERIES AND RECHARGING	11
4.0. MAINTENANCE AND TROUBLESHOOTING	13
5.0. TECHNICAL DATA	14

CONGRATULATIONS ON YOUR VELA BLUES 100

Your new VELA Blues 100 powerchair is designed to make your life easier. We have developed an extremely compact powerchair for indoor use. VELA Blues 100 is designed for simple and logical operation indoors.

When you move about indoors, you will appreciate being able to raise the seat to reach across a table or take something down from a shelf. The seat lift is operated from the powerchair's control box or – depending on the individual construction of the chair – from a separate tilt switch.

We hope that you are happy with your VELA Blues 100. By taking a few minutes to read this manual, you will quickly gain an overview of the construction and possibilities of the VELA Blues 100.

IMPORTANT!

This user manual covers useful information on set-up, operation, maintenance, technical data etc. Please read this manual carefully and keep it for later reference. It must always be available to the user. This user manual must always be kept on the chair.

The manual is also available on VELA's website: www.vela.eu

If you have any questions, please contact our customer service department: VELA, Goeteborgvej 12, 9200 Aalborg, Denmark. Tel. +45 9634 7600.

SAFETY

This chair carries the CE mark, and complies with all current EU regulations.

RECYCLING

Some of the materials in this chair can be recycled, so please return the chair to VELA or to your local recycling centre. Batteries and electronic parts must be returned to VELA or to your local recycling centre.

GUARANTEE AND SERVICE

Your VELA Blues 100 is covered by VELA's 24-month right of return for any problems that can be attributed to material or fabrication defects. The guarantee does not cover parts that suffer wear, such as tyres, batteries, seat cover etc. VELA provides service on weekdays.

Enjoy your VELA Blues 100.

Vermund Larsen A/S

VELA

Goeteborgvej 12

9200 Aalborg

Denmark

VELA reserves the right make changes to this document.

1.0. OPERATION

This powerchair is a class A (EN 12184) electric wheelchair. The first attempts to operate the chair should be made in an area with plenty of floor space.

Choose a low operating speed (press the "down button" a few times) until you have become accustomed to the chair.

JOYSTICK



1.1. BEFORE OPERATION

Before starting the powerchair, you should check the following:

- The footrest is in place and both feet are placed properly on the footplate.
- The electric lift (seat lift) is in the lowest position.
- The electric tilt (seat tilt) is in neutral position.
- Both motors are connected (both connection arms are horizontal).
- The powerchair is turned on by pressing the I/O button.



ENTERING AND EXITING

Always turn off the power before entering or exiting the chair (the I/O button).

FOOTRESTS

Any footrests should be swivelled to the side or the footplate should be folded up before exiting the chair.

Warning! Do not stand on the footplate or footrests – the chair can tip and you may lose your balance.

1.3. OPERATION



VR2 CONTROL BOX

VELA Blues 100 will normally be equipped with a VR2 control box. The powerchair is turned on and off with the I/O button.

To operate a function, press the left button = electric seat lift. To operate a function, push the joystick forwards (up) or backwards (down).

The 'speed up' button increases operating speed: Press the button repeatedly and you will see the indicator light move up.

The 'speed down' button reduces operating speed: Press the button repeatedly and you will see the indicator light move down.

Two of the powerchair's functions are controlled from the control box (normally lift and tilt). Choose one of the functions by pressing the function button, pushing the joystick forwards to move forwards or backwards to move backwards. (It is possible to change between the two functions by moving the joystick to the side). The control box has a signal horn.

STEERING WITH THE JOYSTICK

Move the joystick in the direction you want the chair to move.

Acceleration and speed are controlled by moving the joystick:

Large movement = fast acceleration and high speed.

Small movement = slow acceleration and low speed.

Warning! Do not subject the joystick and control box to vertical pressure and do not push hard on the joystick as this can damage parts in the box!

CONTROLLING THE POWERCHAIR

Pushing the joystick forwards, activates both wheels synchronously and the chair moves forwards. Pushing the joystick backwards, activates both wheels synchronously and the chair moves in the opposite direction. If the joystick is pushed diagonally in one direction, the speed is distributed between the two wheels asynchronously (and in some cases in separate directions). This will cause the chair to move on a path in the right direction.

NOTE! When changing directions, the chair's swivel wheels should change position: they will swing on the axle of the front fork until the wheels are pointing in the direction you wish to go. When the powerchair is standing still before changing direction in this way, a large amount of motor power is needed to drive the wheels. This power can cause a jolt to the front end (leg supports) or back end (if the chair is front-wheel drive) of the chair. The effects of this jolt can be minimised by increasing power gradually and moderately (small movements of the joystick).

Warning! Even though the electrical system is shielded and complies with standard regulations, there is a risk of disrupting the operation of the chair if it is exposed to a field of strong electromagnetic radiation from other devices. In case of such a disruption, the powerchair must be turned off immediately using the I/O button.

NEGOTIATING OBSTACLES

The VELA Blues 100 can handle a difference in surface level of up to 2.5 cm, which makes it possible to cross e.g. doorsteps and many types of ramps. Always drive as straight as possible when negotiating such obstacles. For abrupt changes in level, it helps to increase the length and speed of your approach.

Warning! The electric seat lift must be lowered to the lowest position before negotiating obstacles as a raised seat reduces stability and can cause the chair to tip over!

NEGOTIATING SLOPING SURFACES

The powerchair has enough traction to negotiate standard household ramps. (Max. gradient with load of 75 kg: 12°). When moving down ramps, the chair must always move perpendicular to the sloping surface and at a slow speed.

Warning! The electric seat lift must always be lowered to the lowest position before negotiating sloping surfaces as a raised seat reduces stability and can cause the chair to tip over!

Warning! The chair must never be driven diagonally or sideways across a sloping surface as this can cause the chair to tip over!

BRAKING

To stop the powerchair, return the joystick to the central position or let go of it.

ELECTRIC BRAKE

The VELA Blues 100 is equipped with an electric brake that is activated automatically moments after the drive motors stop. This brake ensures that the chair is fully locked when entering and exiting, parking on ramps etc.

EMERGENCY BRAKE

The powerchair can be stopped immediately by turning it off: Press the I/O button.



Warning! When the powerchair is turned off, the electric brake is activated immediately. This will cause a very sudden stop. The chair must therefore only in emergencies be turned off while it is still moving. There is a risk of the user falling out of the chair and/or of the chair tipping over!

RELEASE

Should it be necessary to manoeuvre the chair manually, the electric brake and motor can be released by turning the two handles located on top of the motor blocks clockwise to the "FREE WHEEL" position.



Note! It is not possible to operate the Blues 100 powerchair electrically when the release function is in use. Turn the two handles back to "DRIVE" immediately after manoeuvring the chair.

Warning! Do not attempt to manoeuvre the chair manually when the motor is still in gear, when the chair is on a sloping surface, on a landing or in a similar situation where the chair might start rolling or lose control.

STANDING STILL AND PARKING

Always use the I/O button to turn the powerchair off as soon as you no longer want to drive it.

NOTE! When the chair is turned on, power is consumed from the batteries. It is therefore a good idea to turn off the chair when it is parked.

2.0. SAFETY POINTS

For your safety, the following should always be observed:

- The VELA Blues 100 is for indoor use only. It should only be operated on level surfaces in places where the user feels confident.
- Always turn the VELA Blues 100 off before entering or exiting the chair, as moving the joystick inadvertently can cause the chair to move.
- Never stand on the footplates when entering or exiting the chair as this can cause the chair to tip over.
- Always lower the electric seat lift to the lowest position before negotiating obstacles and before moving at high speeds.
- Always lower the electric seat lift to the lowest position before negotiating sloping surfaces.
- Never drive sideways or diagonally on a sloping surface as this can cause the chair to tip over.
- The chair must never be subjected to rain or extreme moisture as moisture can damage the electronics.
- Never cover the battery charging unit during recharging as it produces heat. The battery charging unit must be placed on a firm, well-ventilated surface.
- When moving backwards and when operating the electric functions, make sure that other people are well out of the way.

Warning! The VELA Blues 100 is not approved as a passenger seat for use in cars, buses or trains.



3.0. SET-UP

ASSISTANCE WHEN EXITING THE CHAIR

The electric seat lift can help in connection with entering and exiting the chair.

For instance, the electric seat lift can be placed in the highest position before sliding out over the front edge of the seat with legs as stretched out as possible.

The reverse procedure can be followed when entering the chair.



The self-locking conical guideways in the armrest bars ensure that the armrests cannot be pushed down unintentionally. Thus, it is safe to use the armrests as support when entering and exiting the chair.

SITTING POSITIONS

A correct sitting position is achieved when the chair back is placed with approx. three fingers of space between the front edge of the seat and the hollows of the knees. The chair back can be moved (using tools) forwards and backwards above the seat to adjust for seat depth. The armrests are in the correct position when the user is sitting up straight with elbows in to the body and arms resting comfortably without lifting the shoulder blades.



ARMREST ADJUSTMENT

The height of the armrests can be adjusted by loosening the hand screw in the armrest bar.

The distance between the armrest can be adjusted by loosening the hand screw under the seat.

Remember to tighten the hand screws again.



SEAT LIFT

The VELA Blues 100 has, as standard equipment, a seat lift. It is operated by pressing the bottom left function button.

To raise the seat, press the left function button while pushing the joystick forwards.

To lower the seat, press the left function button while pushing the joystick backwards.



FOOTRESTS

ELECTRIC SWIVEL FOOTRESTS, STANDARD

The VELA Blues 100 has, as standard equipment, electric swivel footrests.

The VELA Blues 100 can also be equipped with a set of adjustable footrests attached to the seat frame.

The footrest can be moved out of the way when entering and exiting the chair by lifting the bar approx. 3 cm up and swinging the footrest around the side at the same time.

The depth of the footrest can be adjusted using the tool for adjusting the sliding bar in the seat frame either forwards or backwards.

The footrest can be adjusted to suit leg length by regulating the length bar in its hollow tube. The foot angle can be adjusted using the tool in the swivel part of the footplate.

FOLDING S-FOOTREST WITH PLATE (OPTIONAL)

The footplate saves space and follows the height of the seat when the electric seat lift function is used. The length of the footrest is regulated by loosening the six-sided bolt on the bar down by the footplate. The desired position for the footplate is selected and the six-sided bolt is tightened again. The s-footplate can be folded up when entering and exiting the chair.



ADJUSTMENT OF CHAIR BACK/SEAT DEPTH

The chair back is released by loosening the fastener on the chair back bar located in the back under the seat. The chair back is placed in the desired position and then tightened again. The chair back bar is equipped with a pull stop to ensure that the chair back bar is correctly anchored in the chair back holder. Destruction of the pull stop and extreme pressure on the chair back bar can cause damage to the chair back holder.

PLACEMENT OF CONTROL BOX PARALLEL FITTINGS


The control box is mounted on either the right or the left side according to customer specifications.

The VELA Blues 100 is delivered with parallel fittings as standard, making it possible to slide the control box along the armrest so that the chair can e.g. come all the way up to a table.



4.0. BATTERIES AND CHARGING

BATTERIES

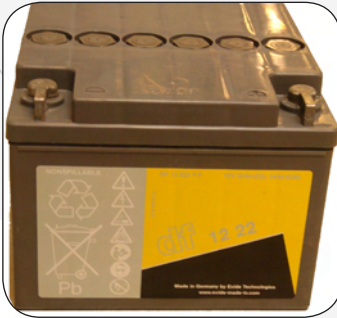

The VELA Blues 100 is equipped with two 12V-24 Ah maintenance-free and vapour-free batteries. Good quality new batteries provide the equivalent of three hours of continuous operation. Experience shows that this covers the daily needs of most users. The batteries are dimensioned to be fully run down 200 times. Their lifetime will be extended significantly if they are only partially run down before charging. Defective batteries must be replaced with batteries of the same type. We recommend that battery replacement be carried out by a service technician.

NOTE! When installing batteries, extreme care must be taken to connect + / - correctly. Red wire to the + pole. Incorrect connection can damage the powerchair's electronics.

NOTE! Batteries must not be disposed of with ordinary household rubbish. Used batteries and other electrical equipment must be deposited at an approved recycling centre according to the marking on the product.




BATTERY CHARGE INDICATOR

The control box has a built in battery indicator. For fully-charged batteries, three groups of coloured LEDs light up. During use, the voltage in the batteries falls and the row of lights decreases: one at a time each light goes off.

- Green signal = DRIVE
- Yellow signal = Recharging needed soon
- Red signal = Recharging needed immediately

The battery indicator will flash to call attention to low battery power.

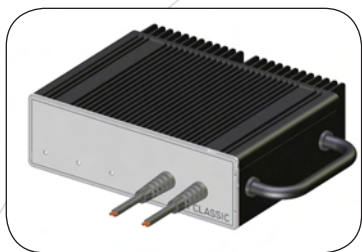

BATTERY CHARGING UNIT

Batteries are charged using the accompanying charging unit (e.g. Classic 8A). The charging time is approx. 8 hours, depending on the condition of the battery. It is recommended that the batteries be charged every night. The charging unit is fully automatic and controls the power and voltage during the charging. When the batteries are fully charged, power consumption will be minimal, so you can safely continue to allow the charging unit to be connected until the next time you use your VELA Blues 100.



WHAT TO DO

1. Park the chair and turn off the power (press the I/O button).
2. Plug the charging unit into the socket on the front of the underside of the control box. **NOTE!** The chair cannot be turned on when the charging unit is connected.
3. Plug the charging unit to the mains power (230 V) and turn it on. The yellow light on the front panel of the charger shows that the batteries are being charged.
4. When the batteries are fully charged, the "completed" light on the front panel of the charger will be green.

NOTE! For extended periods of disuse (days/weeks/months) the battery charger should be connected to ensure optimum power when your VELA Blues 100 is used again.

NOTE! In some cases, batteries that have been run down too far (flat) cannot be recharged by the charging unit. This is for your own safety, because dangerous vapours can develop when flat batteries are recharged. Call your VELA service technician.

4.0. MAINTENANCE AND TROUBLESHOOTING



4.0.1. CLEANING

Wipe/wash the wheels as needed using mild soap and water. Wipe the protective motor cover with a soft cloth. Do not use abrasive sponges or brushes.

Clean the seat and back with a vacuum cleaner or a brush. The fabric is a standard wool type, which has been chosen for optimum comfort. It is not recommended that the seat cover be washed as it may shrink. Dry-cleaning recommended.

The control box can be cleaned with a damp cloth and a little washing-up liquid if necessary. The armrests can be wiped off with a damp cloth and warm soapy water if necessary.

4.0.2. TROUBLESHOOTING

If your VELA Blues 100 should stop working, check the following first:

CHAIR DOES NOT START

- Speed indicator is blinking: Check that both connection handles above the motor are set at "DRIVE".

CHAIR LOSES POWER

- Check that the charging unit is connected correctly.
- Check that the charging functions normally (read LED lights on front panel of charger).
- Call a service technician to have the condition/age of batteries checked.
- Call a service technician to have the wiring to batteries checked.

7.0. TECHNICAL DATA

CHAIR IDENTIFICATION NUMBER

The serial number is found on the chassis between the swivel wheels.

DRIVE MOTORS

Two 24V Leroy motors 1500 rev./min.

Electronic disc brakes.

BJ Gear worm drive.

BATTERIES

Two 12V-24Ah vapour-tight and maintenance-free powerchair batteries.

Dimensions: 125 x 167 x 176 mm

Charging time approx. 8 hours with the battery charging unit.

ELECTRONICS

Power network: 24VDC analogue wiring

PG drivers powerchair electronics VR2

Electric seat lift: Actuator VELA Linak LA30.1.200.24VDC

CHARGER

Automatic Power Charger Classic 8A / 230 VAC

7.0.7. TECHNICAL DATA

	VELA Blues 100
Drive motors	2 x 24V
Locking brake	24V DC Electromagnetic
Batteries	2 x 24 Ah
Electronics	VR2
Charger	220-230 24V/Power charger 8A classic
Drive wheels	200 x 50 mm solid
Swivel wheels, front	175 mm PU
Swivel wheels, back	-

Seat size	44 x 36 cm, 44 x 44 cm, 46 x 40 cm, 49 x 46 cm
Sitting height	50-70 cm
Seat height front edge, no-load	50 cm
Height without seat	40 cm
Back size	31 x 16 cm, 38 x 22 cm 38 x 26 cm, 39 x 29 cm, 38 x 40 cm
Back tilt, manual	20°
Back tilt, electric	20° (95°-115°)
Clearance	4 cm
Frame, width	55 cm
Max. load	110 kg
Weight	66 kg

Max. speed	4.7 km/h
Range	22 km
Max. gradient	12°
Obstacles, max.	3.5 cm
Turning radius	75 cm
HMI no.	10900 rear-wheel drive 24936 front-wheel drive

